




UroBilanz User Manual

A clear guide to your local urine and fluid record

UroBilanz helps you organize and review local record data. This manual explains every visible area, the important controls, and a useful first workflow.

English | July 2026 edition

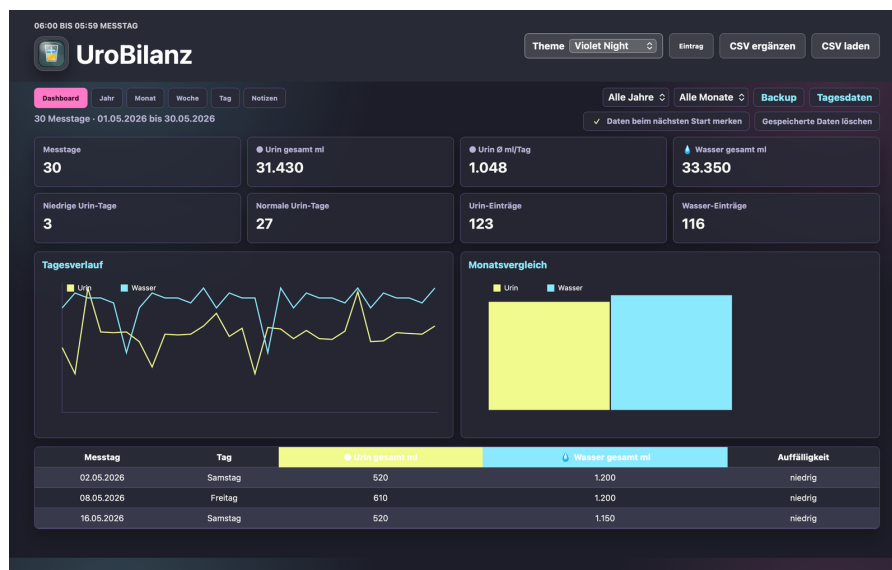

UroBilanz - web app with synthetic demo data


Contents

1. UroBilanz at a glance
2. Where to find each function
3. First useful start
4. Add and manage entries
5. Views and evaluation rules
6. Report, backups and exports
7. Themes, language and settings
8. Privacy, AI help and bug reports

1. UroBilanz at a glance

UroBilanz is a local logging and review tool for urine, water, and note entries. It does not create diagnoses or provide medical recommendations. The web app and macOS app use the same core functions.


- Load CSV: import an original Urinote CSV or daily-data CSV locally.
- Entry: record urine, water, or a note manually.
- Views: compare Dashboard, Year, Month, Week, Day, and Rules.
- Keep copies: export data locally as a backup, daily data, or JSON.

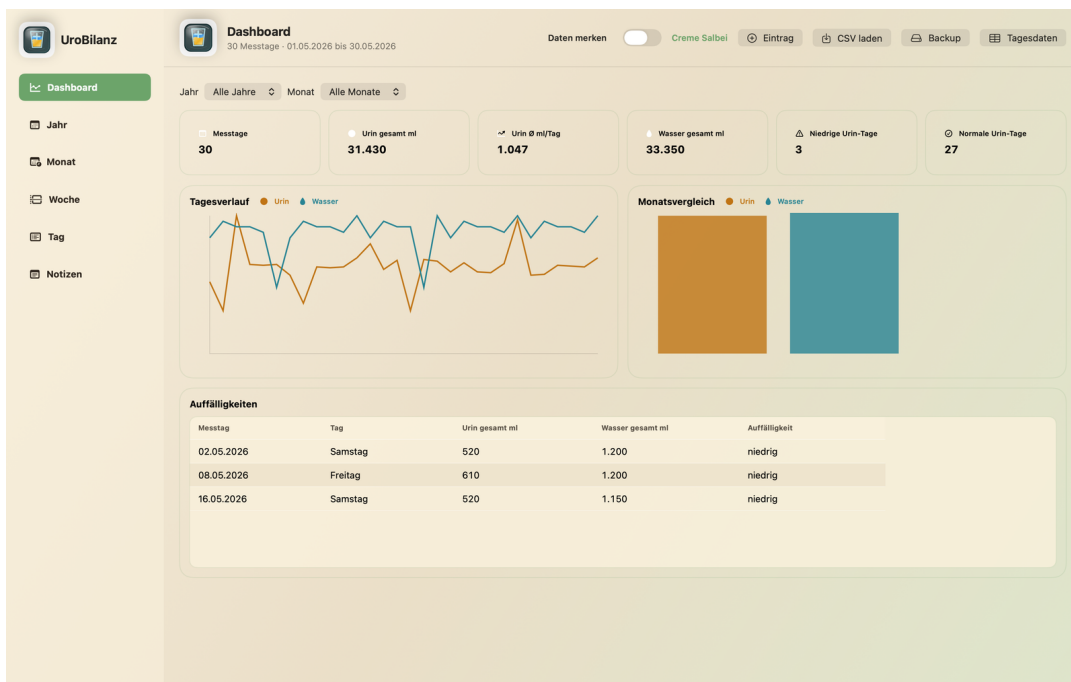

SwiftUI app with demo data

2. Where to find each function

Function	Location and effect
Load CSV	Top right and in the empty start view. Opens a local file dialog; the file stays on your device.
Merge CSV	Top right. Adds new entries from another original Urinote CSV and skips entries already present.
Entry	Top right. Opens manual entry for urine, water, and notes.
Dashboard / Year / Month / Week / Day / Rules	Tabs on the left of the toolbar. They only change the local presentation of loaded data.
Year and Month	Filters on the right of the toolbar. Limit visible tables and charts to the selected period.
Backup	Menu on the right: save a complete backup, daily backup, or JSON export locally.
Medical report	Right side of the toolbar. Choose a period and contents for a local print-friendly report.
Remember data / Delete saved data	Below the toolbar. Saves the loaded CSV only locally in the browser or macOS app, or removes it again.
Theme and Language	Top right. Change appearance or visible language, never measurement data.

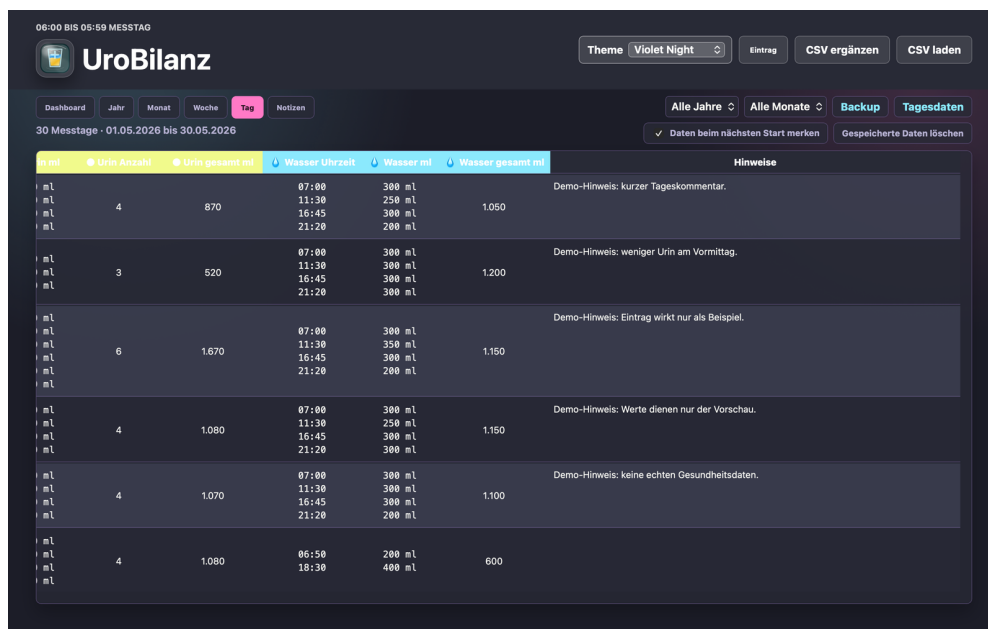

Tabs, filters, and day view

3. First useful start

The app deliberately starts empty. This keeps existing data from other people or app variants separate.

- 1. Open UroBilanz. As long as no data is loaded, you see the start view.
- 2. Choose Load CSV and select an original Urinote CSV or a daily-data CSV.
- 3. Review the Dashboard and switch to Day or Week when useful.
- 4. Enable Remember data only when the CSV should stay stored locally on this device.
- 5. Save important states regularly through Backup.

Merge CSV is only for another original Urinote CSV. A daily-data CSV is an export format and cannot be merged.


The screenshot shows the UroBilanz web application interface. A modal window titled "Eintrag hinzufügen" (Add entry) is open, allowing manual entry of data. The form includes fields for "Datum" (Date), "Urin Uhrzeit" (Urine time), "Urin ml" (Urine volume), and "Wasser Uhrzeit" (Water time). Below these, there are sections for "Wasser ml" (Water volume) and "Tagesnotiz / Hinweis" (Daily note / Note). The background shows a dashboard with a table of entries for May 2026, including columns for "Jahr", "Monat", "KW", "Messtag", "Urin", and "Wasser". The interface is dark-themed with a "Violet Night" theme selected.

Manual entry in the web app

4. Add and manage entries

Entry opens manual input. Date and time determine which measurement day receives a value. Enter at least a urine value, water value, or note. Add keeps the dialog open; Add & close closes it afterwards. Reset clears the dialog. In the Day view, you can edit or delete individual entries; Delete day removes all entries of that measurement day from local analysis and future exports.

- Urine ml and Water ml: amounts stay separate; water is not added to urine volume.
- Note: visible context, but no automatic medical assessment.
- Urine note: a note for a urine value stays aligned with that time.

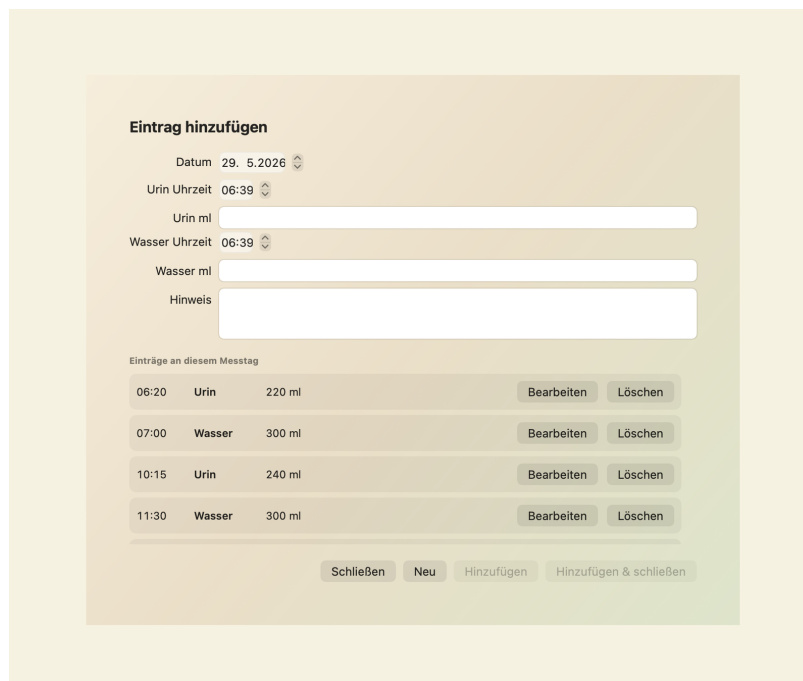

Eintrag hinzufügen

Datum 29. 5.2026

Urin Uhrzeit 06:39

Urin ml

Wasser Uhrzeit 06:39

Wasser ml

Hinweis

Einträge an diesem Messtag

06:20	Urin	220 ml	Bearbeiten	Löschen
07:00	Wasser	300 ml	Bearbeiten	Löschen
10:15	Urin	240 ml	Bearbeiten	Löschen
11:30	Wasser	300 ml	Bearbeiten	Löschen

Schließen Neu Hinzufügen Hinzufügen & schließen

Manual entry in the macOS app

5. Views and evaluation rules

The Dashboard summarizes the loaded period. Year, Month, and Week condense the same data into tables. Day shows individual times, amounts, and notes. Rules explains the organizational logic in use. Streak counts consecutive measurement days with entries.

Rule	Meaning
Measurement day	Runs from 06:00 to 05:59. Entries from 00:00 to 05:59 belong to the previous day; the time remains visible.
Complete	Only days with at least eight hours between the first and last urine or water value are included in totals, averages, and evaluation.
Low	A complete day with less than 700 ml of urine is shown as low.
Normal	A complete day with 700 ml or more of urine is shown as normal.
Incomplete	Boundary days stay visible but are excluded from totals, averages, and low/normal evaluation.
Water and notes	Stay visible separately as context and do not change evaluation.

The low, normal, and incomplete labels are organizational evaluations. They are not a diagnosis and do not replace medical advice.

6. Report, backups and exports

Open Backup in the toolbar and choose the suitable format. Every file is saved at a location you choose through the browser or macOS save dialog.

Option	What it creates locally
Complete backup	CSV with all entries in the original format for later restoration or review.
Daily backup	CSV with summarized measurement days for exchange or archive use.
JSON export	Structured export of entries for your own local processing.
Medical report	Print-friendly report with period, summary, daily progress, and optional daily details and notes.

For the medical report, choose the period first. Daily details and notes are separate switches: leave them off when they are not needed for the intended purpose.

Share exported files only deliberately. An export can contain the content you recorded.

7. Themes, language and settings

The Theme menu at top right switches between built-in themes. Import theme reads a local JSON theme file. Export theme creates an editable local copy of the current theme. Delete theme is active only for imported themes; built-in themes remain available. The DE/EN language choice changes all visible text, including first-start help. In the macOS app, Settings also contains the update check; it is started deliberately and checks only public release information.

- Remember data: stores the currently loaded CSV locally for the next start.
- Delete saved data: removes this local remembered copy; the original file outside the app remains unchanged.
- About UroBilanz: shows version, license, project page, and contact.

8. Privacy, AI help and bug reports

UroBilanz processes measurement data locally. There is no cloud sync and no data backend. Visible GitHub, Discord, license, and contact links open only after a click. Report bug creates a reviewable text draft; CSV values, notes, and health data are not added automatically.

In the empty start view, you can open the manual or choose ChatGPT, Google Gemini, or Claude. Only your click copies a fixed general question with this public manual link to the clipboard and opens the selected service. The app does not pass entries, measurements, file paths, or saved data. Paste the question yourself with Cmd+V and send it only when you choose to do so.

For symptoms, uncertainty, or changes in your health, seek medical advice. This manual and the app do not replace professional guidance.